


## Fluorospot H

# RX

Copyright © Siemens AG 1994

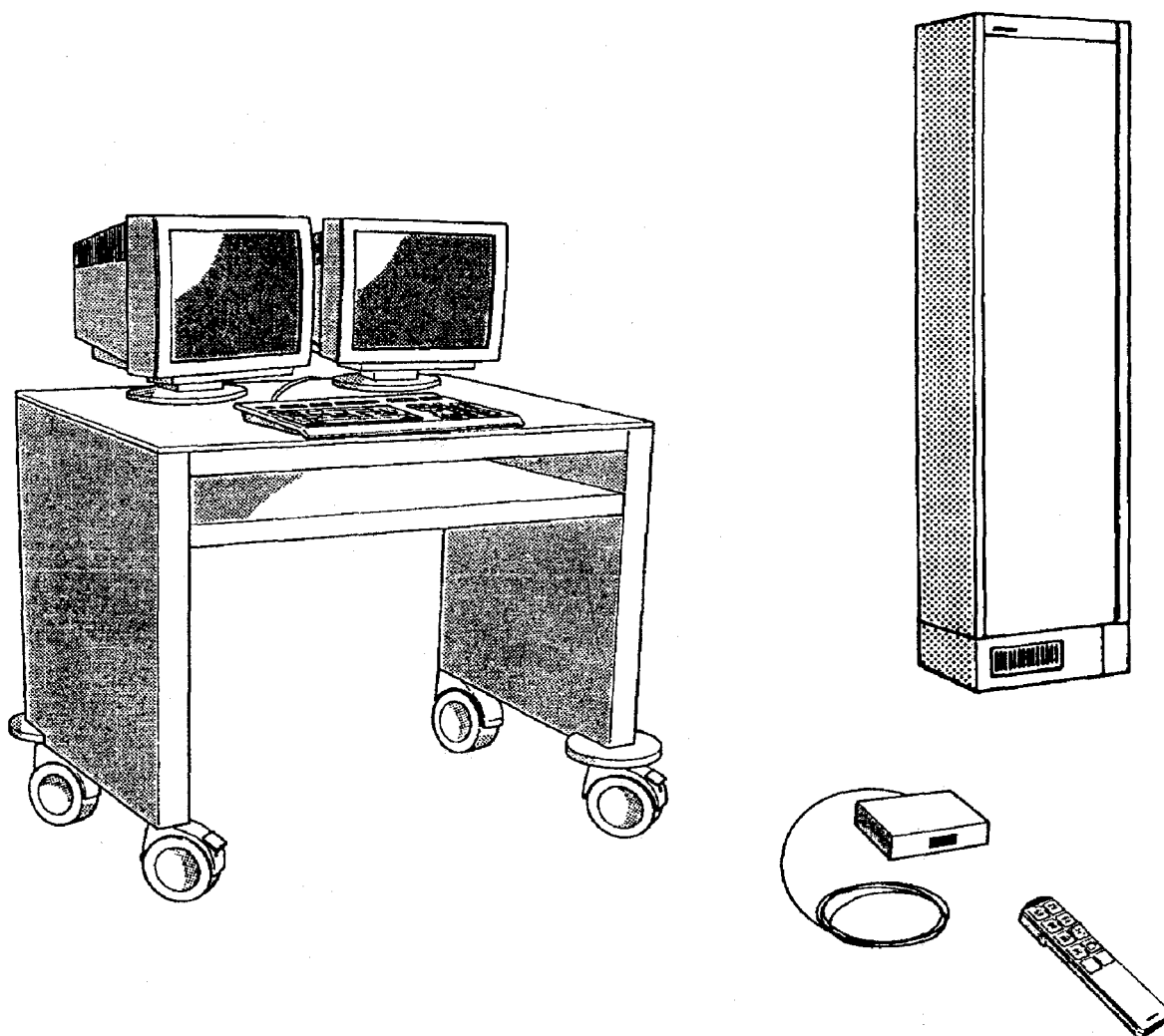

Planning guide


### Register 3

RX41-020.021.01.04.02

Replaces: RX41-020.021.01.03.02

11.95

English

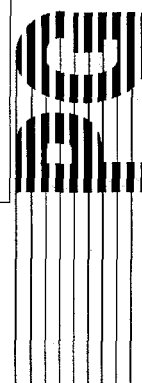

Weitergabe sowie Vervielfältigung dieser Unterlage, Verwertung und Mitteilung ihres Inhalts nicht gestattet, soweit nicht ausdrücklich zugestanden. Zuwiderhandlungen verpflichten zu Schadenersatz. Alle Rechte für den Fall der Patenterteilung oder GM-Eintragung vorbehalten.

Proprietary data, company, confidential. All rights reserved.  
Confié à titre de secret d'entreprise. Tous droits réservés.  
Confiado como secreto industrial. Nos reservamos todos los derechos.



	Page	Rev.
Cover page	0-1	04
List of contents, Revision status	0-2	04
<b>1 Planning</b>		
Planning notes	1-1	03
Dimensions Electronics cabinet	1-2	03
Dimensions Control and evaluation console	1-3	03
Dimensions IR remote control	1-4	03
<b>2 Preparation for installation</b>		
On-site electrical installation	2-1	04
<b>3 Technical Data</b>	3-1	04
<b>4 Electrical Installation</b>		
Notes on cable routing	4-1	03
Protective ground wire routing	4-2	04
RAS wiring	4-3	03
Overview of fixed points	4-4	03
List of fixed points	4-5	04
List of fixed points	4-6	03
<b>5 Changes as compared to previous version</b>		
Changes as compared to previous version	5-1	01

### General notes

**With the distribution of revision level 04, all preceding Speed Infos and drafts are invalid.**

The provisions of the relevant fire protection regulations must be observed for the premises.

When carrying out planning and project organization, the data found in file R A0 must be taken into consideration.

All layouts issued by the Planning Departments must bear a note referring to the installation and delivery conditions of Siemens Medical Engineering Group. The installation and delivery conditions must be submitted with the layouts.

Löchel	TD RX 1	Tel. 09131/84 - 5428	
Weller	TD RX 4	Tel. 09131/84 - 8239	Hotline 7774

### General notes

- All power lines must be run with shielding.
- Shielded power lines are urgently recommended.
- Run protective ground wire parallel with power lines.
- The shielding of the power line must be connected in the electronics cabinet over the shortest possible path to potential ground (PE).

### Power lines belonging to other systems

- Must be run separate (at least 100 cm distance).
- Shielded leads must be used in each case.



### Power line connection and PE cabling

To prevent possible interferences


- Observe the building power line connection, see page 2-1.
- Observe the protective ground wire cabling, see page 4-2.

### Note:

Fluorospot H and additional equipment are not designed for use in areas where there is an explosion hazard.

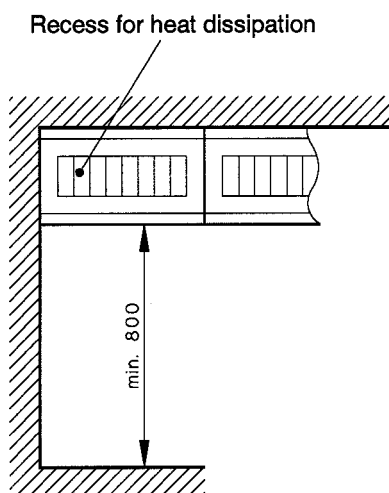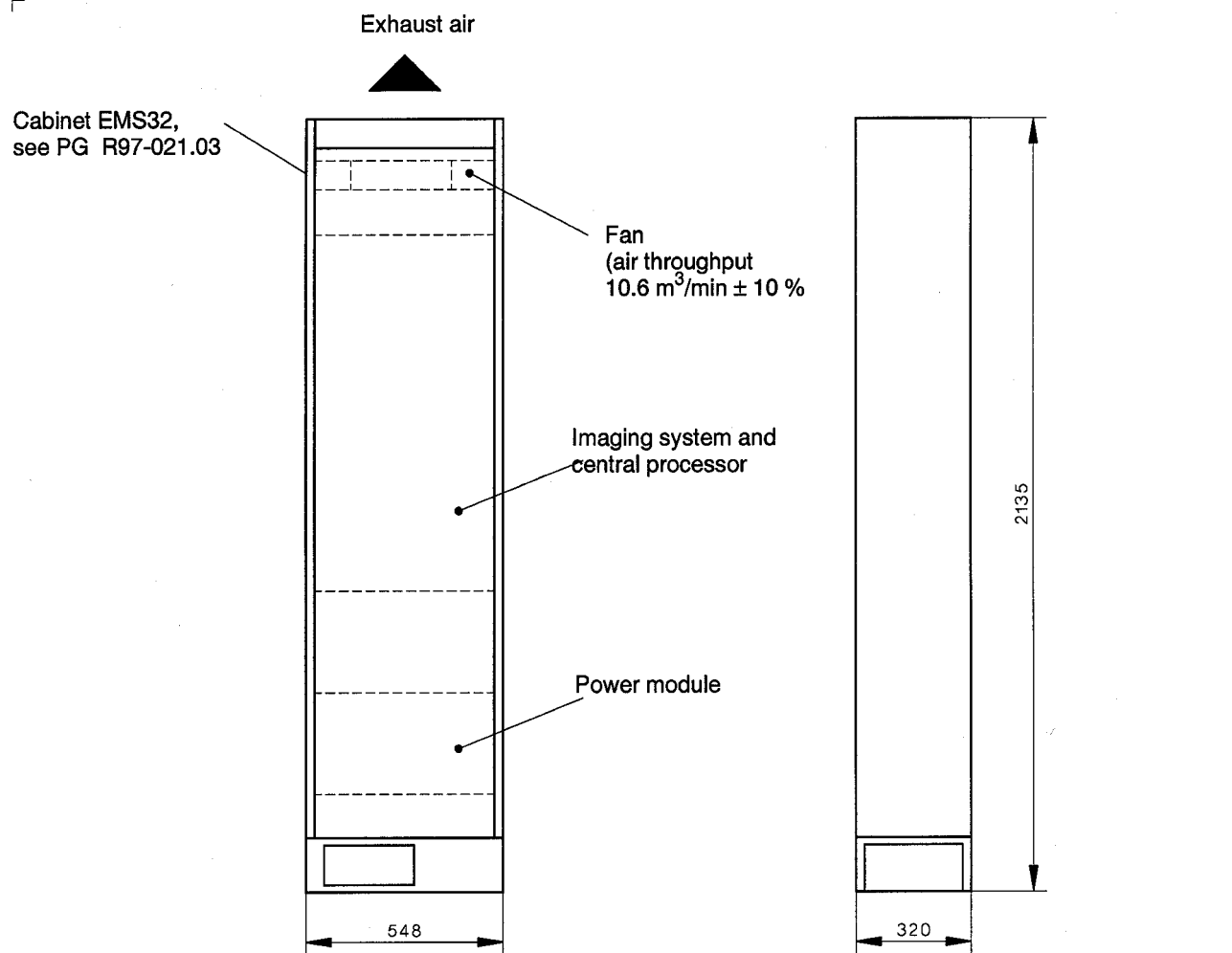


Floor fastening

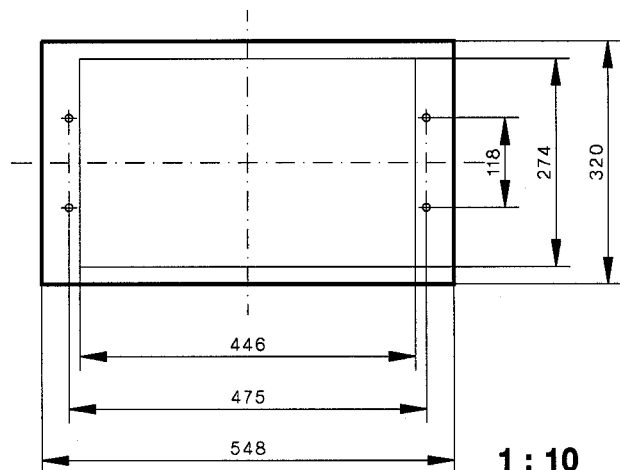**Note :**


The Fluorospot H electronics cabinet can be installed together with further EMS 32 cabinets. However, 1 m distance must be kept to the generator power cabinet at least.

e. g. with monitor table (approx. 45 kg without monitor)

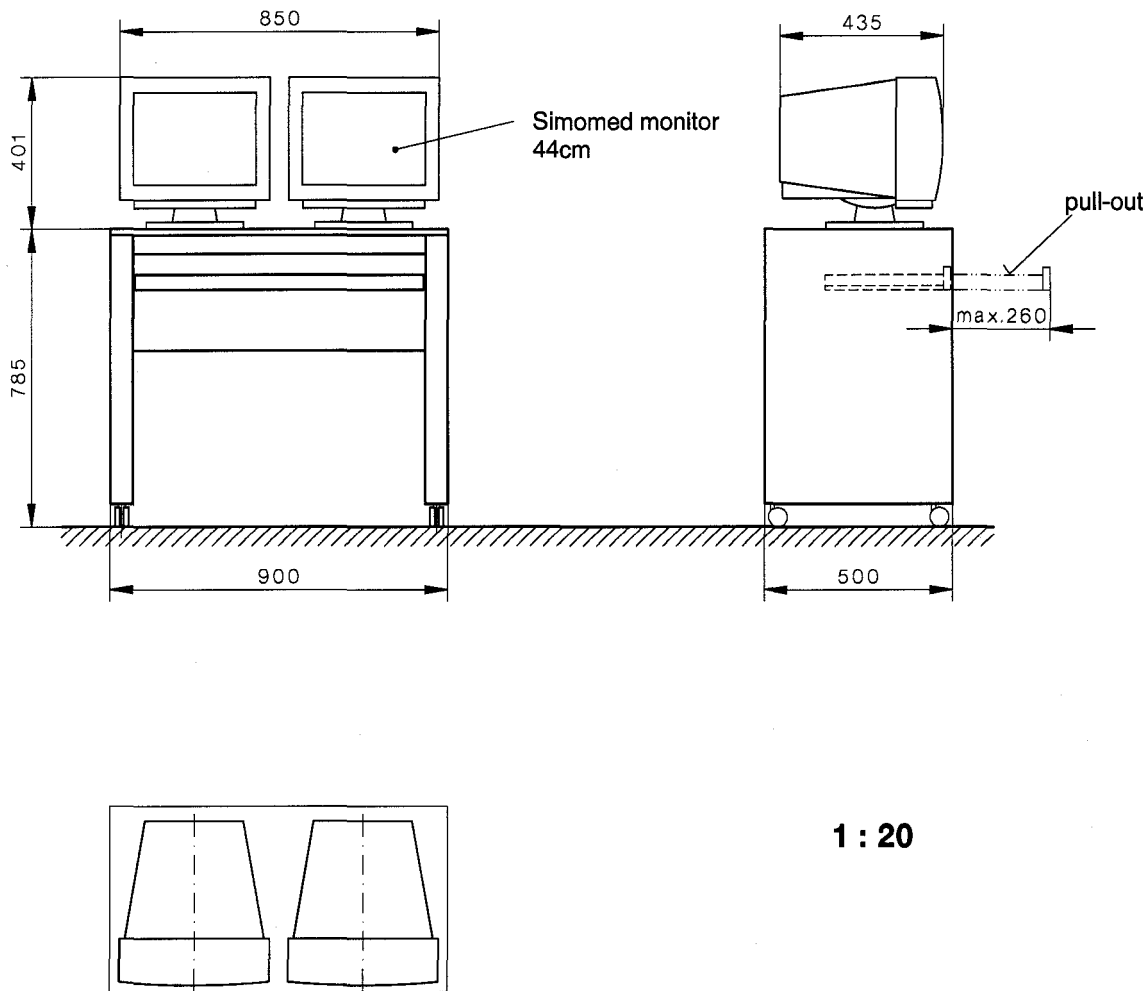

### Note:


The mobile monitor table is an option.  
The customer can choose his own table.

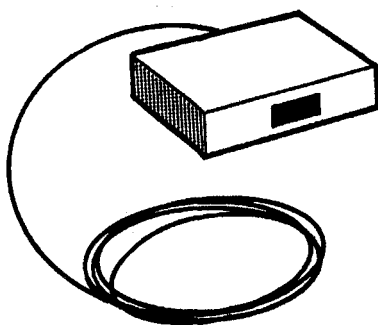

25 m or 50 m lang  
included in delivery

### IR - receiver and preamplifier

Casing:	IP 65, Color: light grey
Supply voltage:	12 V DC
Current consumption:	10 mA
Sensitivity with all Volltronic 3500-transmitters:	50 m when the transmitter is aligned to $\pm 10^\circ$ ambient light $\leq 1000$ lux
Range reduction at sun light:	max. 50 %

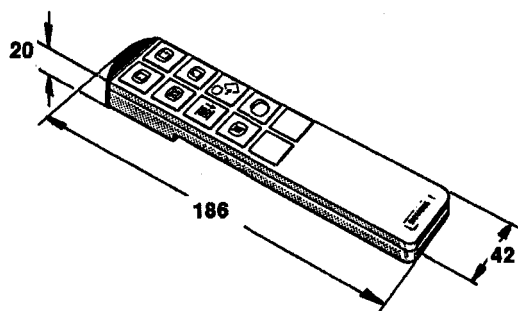


### IR - hand-held transmitter

Casing:	Color: black 42 x 186 x 19 mm
Range:	approx. 50 m with direct emission of receiver $\pm 5^\circ$ and ambient light $\leq 1000$ lux

Charing socket for 12 V - battery, can be charged with 18 V charger

### Proposal for customer's power connection

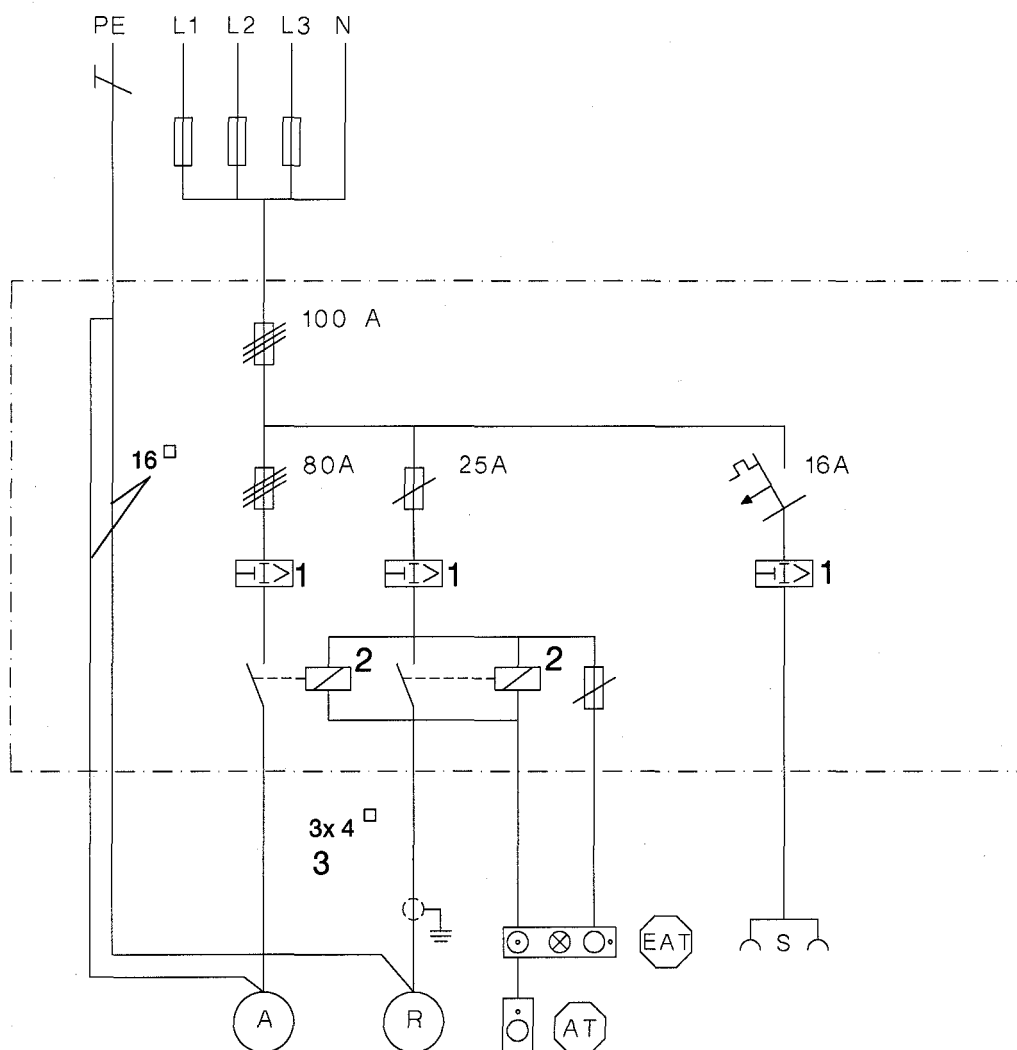

- (R) Fluorospot H Electronics cabinet
- (A) Generator
- (AT) Emergency off button with interlocking circuit
- (EAT) On-Off button with pilot lamp
- (S) grounded socket for Service

- 1 Fault current protective switch  $I_{\Delta N}$  30mA for alternating and pulsating direct currents (surge current proof)
- 2 System contactor
- 3 Note! Connect R via shielded power cable only.

Electrical data	Control cabinet	monitor 44 cm
Power line connection	1 ~ 230V                    + 10/-15 % 50 Hz                        ± 5 %  1 ~ 208, 254, 277V    + 10/-15 % 60 Hz                        ± 5 %	1 ~ 115, 230V            + 10/-15 % 50/60 Hz                   ± 1 Hz
Fuse internal	25 A	2 A flink
Power consumption	max. 1.5 kVA	approx. 140 VA
Heat dissipation	max. 1.5 kW	approx. 90 W

Ambient conditions acc. to DIN IEC 601-1	Control cabinet	monitor 44 cm
perm. ambient temperature during operation during storage	+10° - +40° C - 40° - +70° C	+5° - +40° C - 25° - +52° C
perm. relative humidity during operation during storage	15 - 80 % 10 - 90 %	15 - 75 % 10 - 75 %

Weights and noise development	Control cabinet	monitor 44 cm
Noise development	max. 57 dBA	-
Weight	max. 186 kg	approx. 22 kg

Surface color	
Main color	white mottled lacquer; Med surface No. 4146 similar to RAL gray-white 9002
Combination color	plain gray mottled lacquer, Med surface No. 4426 similar to RAL dust-gray 7037

Packing	
largest crate L x B x H	2380 mm x 880 mm x 920 mm
Weight	286 kg

### Notes on cable routing

Minimum depth of cable duct shall be 60 mm  
In case of cable crossings, a larger duct might be necessary.


Separate the high voltage cable and power supply leads from the control cables and the video signal leads.  
(if possible, shielding measures should be provided).

The cable list shown below provides a subdivision into four categories in accordance with technical safety and functional point of view.

Category 1 (more than 300V)	Category 2 (max. 300 V)	Category 3	Category 4
High voltage cable	Power cables	EK14	Signal cable
Rotating-anode cable	EK20	EK48	Video cable
Protective ground wire	EK34		Fiber optic cable
Water hose			


### Proposal for cable routing


- Cable routing in separate conduits or closed cable ducts
  - The minimum diameter results from the cable list, column "minimum passage"
  - The cable with the largest plug must be drawn in first
- Separate cable routing in open cable ducts with metal dividers or similar

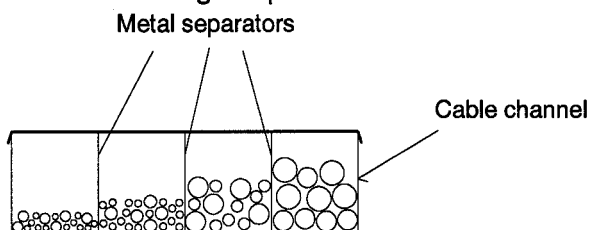

- Calculation of the minimum cross-section  
 $\Sigma$  of the existing cable cross-sections according to the cable list, column "**required cable duct**"
- Specified bulk factors have been taken into account


### Power distribution box

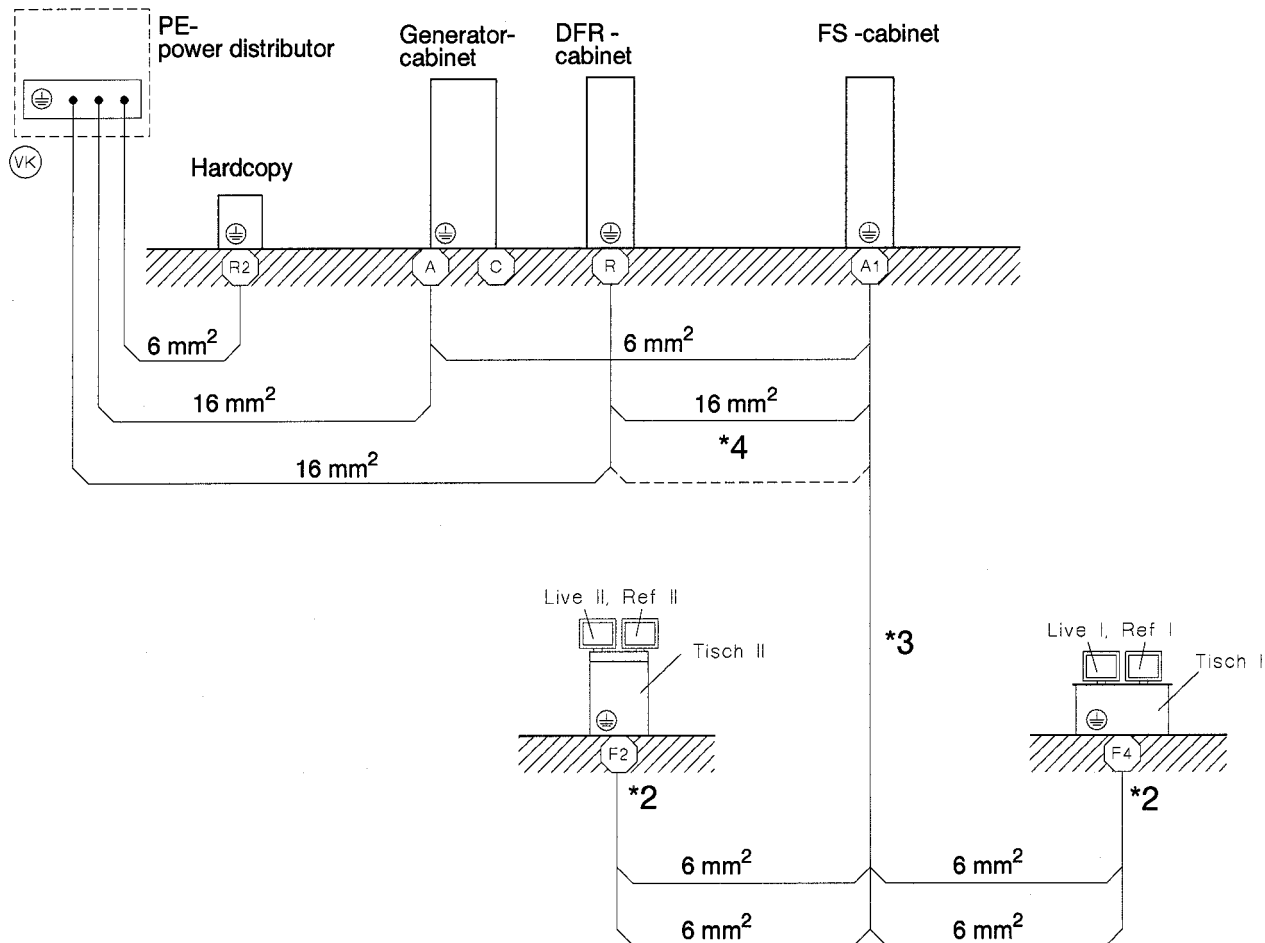


- \*1 See installation and setting instructions VID - H1X R52-040.033 Page 4-3
- \*2 With monitors X1763/X1765, disconnect the jumper between video ground and 
- \*3 Not applicable for Simomed monitors. With Simomed monitors, grounding via power line cable
- \*4 Remove connection, if present.

### RAS - connection (LWL)

Connect fiber optic cables as shown below.

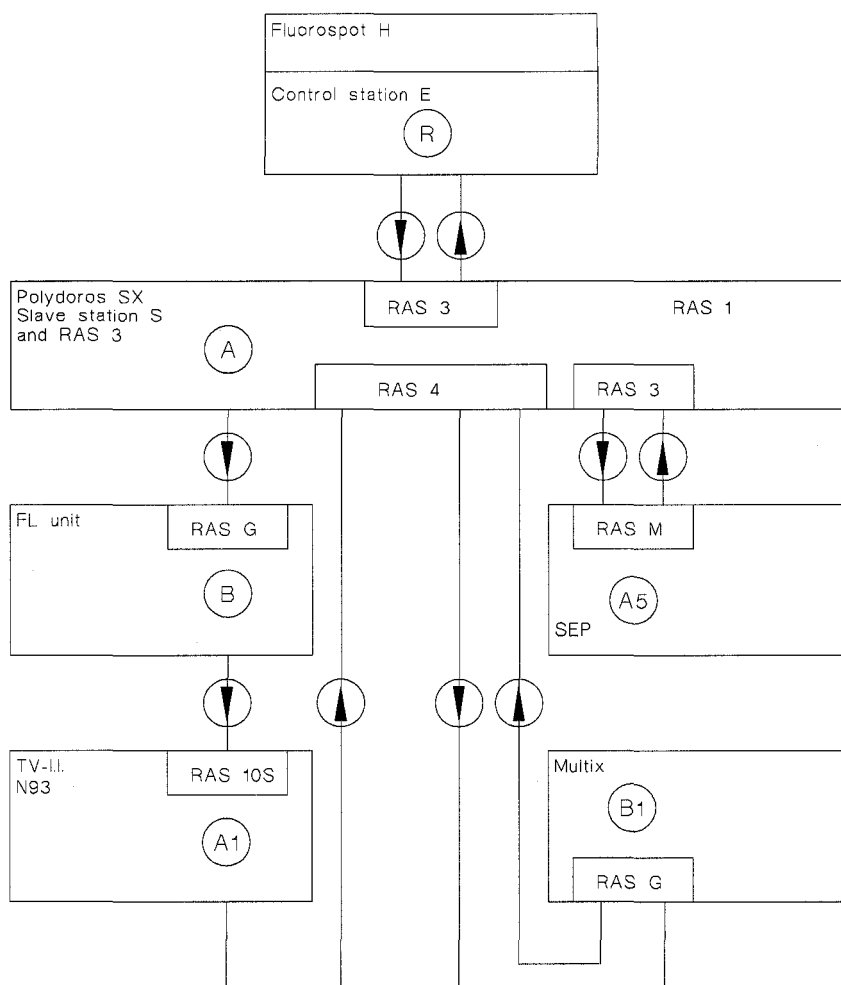

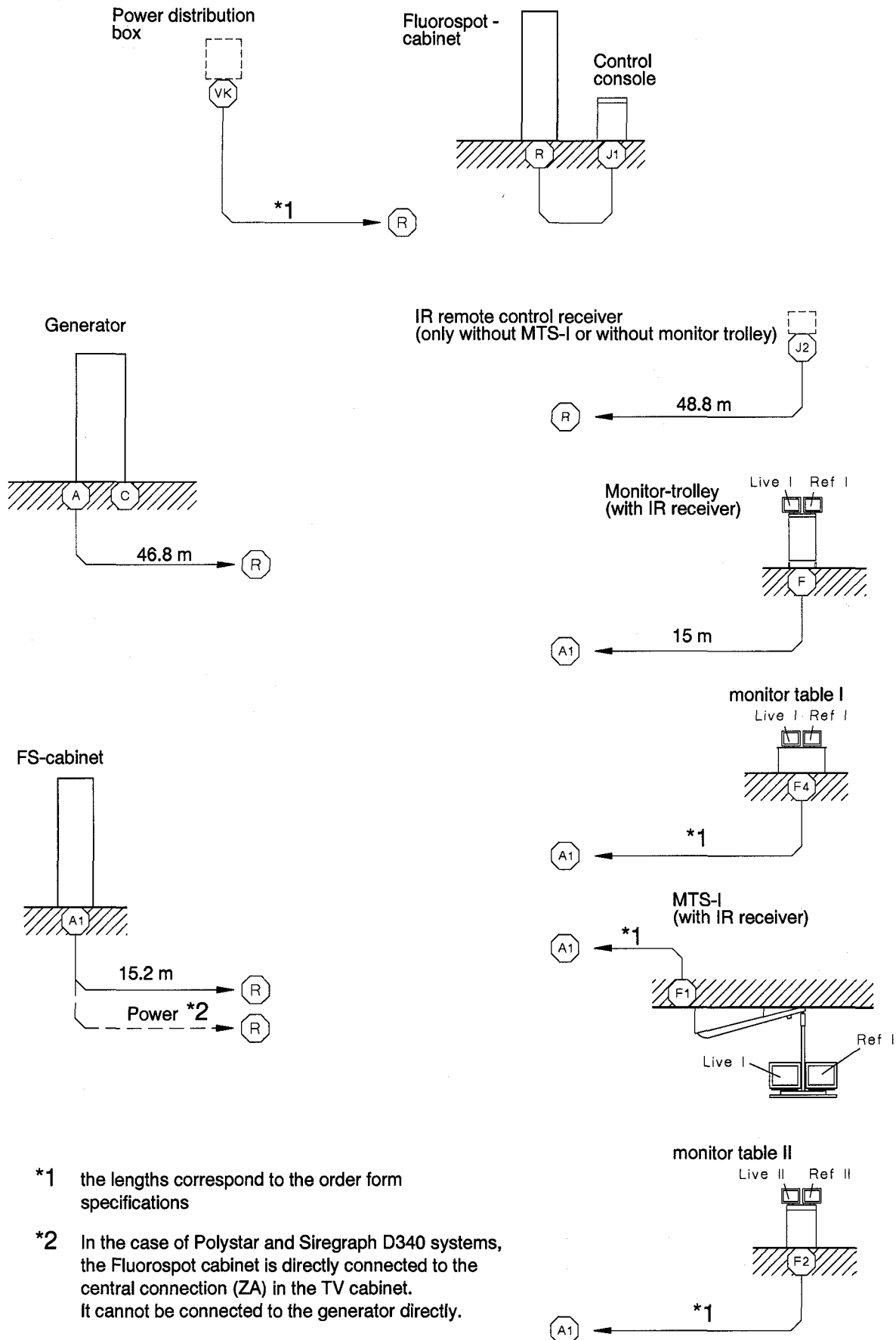

## List of fixed points

from fixed point	to fixed point	Cable type and category ( )	cable dia [mm]	min. opening Ø [mm]	required cable duct [mm <sup>2</sup> ]	Remarks
---------------------	-------------------	--------------------------------	-------------------	------------------------	---	---------

### Cabling FL-H

internal cable length[m]	2.0	A	R	1.2	internal cable length[m]						
2.0	A	R	1.2			LWLG	(4)	8.0	11	80	RAS
2.0	A	R	1.2			LWLG	(4)	8.0	11	80	RAS
0.8	R	VK	0.0			NK 3x4 mm <sup>2</sup> a(1)	19	20	451	Power cable	
0.8	R	VK	0.0			PE16 mm <sup>2</sup>	(1)	19	15	97	
1.6	A1	R	0.5			PE16 mm <sup>2</sup>	(1)	8.8	15	97	
1.2	A1	R	1.2			SK111	(4)	8.0	60	80	
1.6	A1	R	1.2			TriBH	(4)	7.5	30	71	
1.2	R	J1	0.0			SPK	(4)	7.0	40	62	*1

### Cabling IR remote control receiver

1.2	R	J2	0.0			IRSpe	(2)	10	40	125	only with 9816
-----	---	----	-----	--	--	-------	-----	----	----	-----	-------------------

### Cabling Monitor-trolley

1.5	A1	F	1.5			SGK	(2)	22	25	605	monitor trolley
-----	----	---	-----	--	--	-----	-----	----	----	-----	-----------------

### Cabling Monitor-trolley table I

1.5	A1	F4	1.5			NK 3x0.75mm <sup>2</sup>	(2)	8.5	60	91	monitor table I
1.5	A1	F4	1.5			NK 3x0.75mm <sup>2</sup>	(2)	8.5	60	91	monitor table I
1.5	A1	F4	1.5			PE6 mm <sup>2</sup>	(1)	7.3	15	67	monitor table I
1.5	A1	F4	1.5			PE6 mm <sup>2</sup>	(1)	7.3	15	67	monitor table I
1.5	A1	F4	1.5			Triax	(4)	7.5	50	70	monitor table I
1.5	A1	F4	1.5			Triax	(4)	7.5	50	70	monitor table I

\*1 supplied with along Fluorospot

a = shielded

from fixed point	to fixed point	Cable type and category ( )	cable dia.[mm]	min. opening Ø [mm]	required cable duct [mm <sup>2</sup> ]	Remarks
---------------------	-------------------	--------------------------------	-------------------	------------------------	---	---------

**Cabling Monitor-trolley table II**

internal cable length[m]	1.5	A1	1.5	F2	internal cable length[m]	1.5	NK 3x0.75 (4)	8.5	60	91	monitor table II
1.5	A1	1.5	F2	1.5	NK 3x0.75 (4)	8.5	60	91	monitor table II		
1.5	A1	1.5	F2	1.5	PE 6 mm <sup>2</sup> (1)	7.3	15	67	monitor table II		
1.5	A1	1.5	F2	1.5	PE 6 mm <sup>2</sup> (1)	7.3	15	67	monitor table II		
1.5	A1	1.5	F2	1.5	Triax (1)	7.5	50	71	monitor table II		
1.5	A1	1.5	F2	1.5	Triax (1)	7.5	50	71	monitor table II		

**Cabling MTS-I**

1.5	A1	1.1	F1	NK 2x1.5mm <sup>2</sup> (2)	9.0	20	102	MTS-I
1.5	A1	1.1	F1	PE 16 mm <sup>2</sup> (1)	8.8	15	97	MTS-I
1.5	A1	1.1	F1	SPK (2)	14.5	70	263	MTS-I
1.5	A1	1.1	F1	Triax (4)	7.5	50	71	MTS-I
1.5	A1	1.1	F1	Triax (4)	7.5	50	71	MTS-I

Page	Changes
------	---------

0-1	extended
0-2	extended
2-1	FI description changed
3-1	Technical data adapted
4-2	text changes
4-5	cabling changed
5-1	newly inserted
5-2	newly inserted

This page has been left blank intentionally.